

CERTIFICATE

Date of Manufacture _____

Inspection Result *qualified*
QC002

Inspector No _____

INTENSE PULSED LIGHT

HAIR REMOVAL DEVICE
USER GUIDE




TABLE OF CONTENTS

Device description.....	1	Before use.....	8
How (IPL) Intense Pulsed Light Device works ?..	2	1.Get your skin ready.....	8
		2.Skin test.....	8
		3.Hair test.....	8
Contraindication.....	3	While use.....	9
1.General conditions.....	3	1.Selecting the right intensity.....	9
2.Medications/History.....	4	2.Handling the device.....	10
3.Disorders.....	4	3.Treatment schedule	11
4.Skin condition.....	5		
5.Areas.....	5		
		After use.....	12
Important notes.....	6	1.After care.....	12
1.Warning.....	6	2.Common reactions.....	12
2.Device Care.....	6	3.Rare side effects.....	12
3.Caution.....	6		
		Device cleaning and storage.....	13
Course of treatment.....	7	Device Specifications.....	13
1st Treatment - (Week - 0).....	7	Frequently asked question.....	14
2-3 Treatment - (Week - 2 & Week 4).....	7	Device Warranty Card.....	16
4-5 Treatment - (Week - 6 & Week 8).....	7		


Thank you for purchasing this product. In order to facilitate the use, please read this manual carefully and keep it properly. The pictures in this manual are for illustration only.

Device Description


1.1. Instrument description

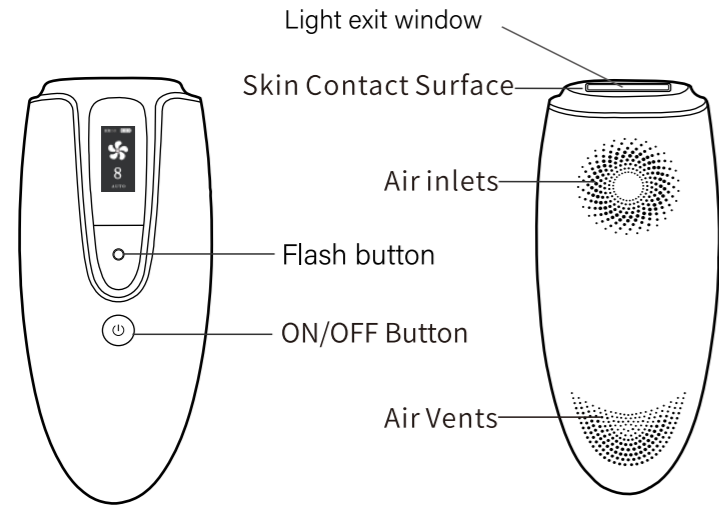

1.2. IPS display interface

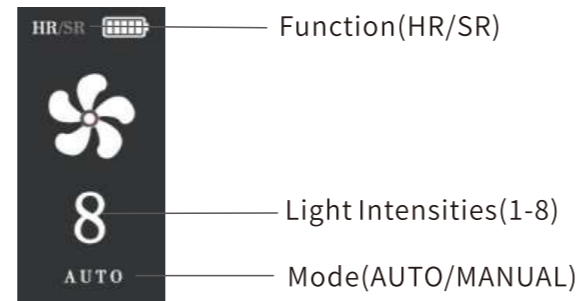


1.3. Random accessories

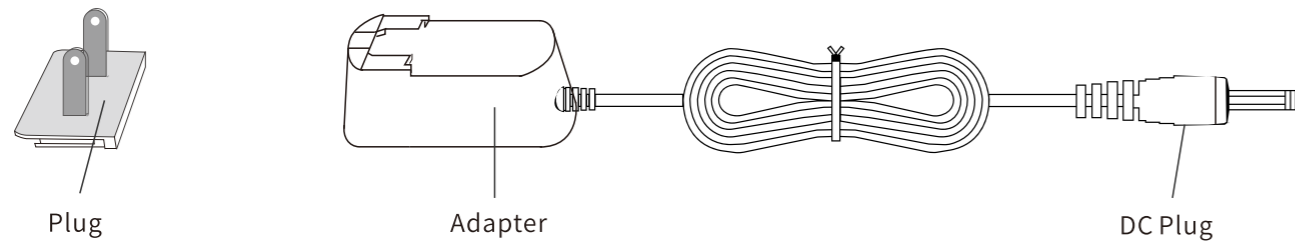


How (IPL) Intense Pulsed Light Device Works


The IPL (Intense Pulsed Light) device projects light pulses (with wavelength range between 550nm-1200nm) on to the skin. During the process melanin in the hair absorbs the light pulses emitted from the device, thereby heating and destroying the hair cells. This painless process then stimulates the hairs to go into the resting phase. The device achieves long-lasting hair-free skin.

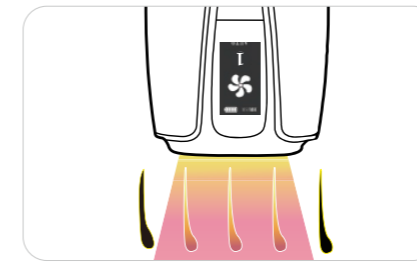


The darker the hair, the better the pulses of light are absorbed.

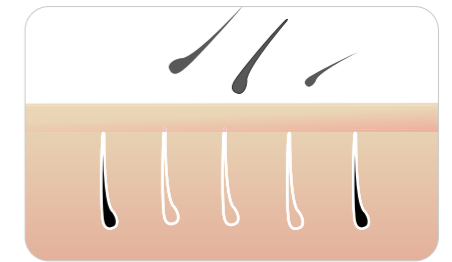

The pulses of light stimulate the hair follicle to go into a resting phase.

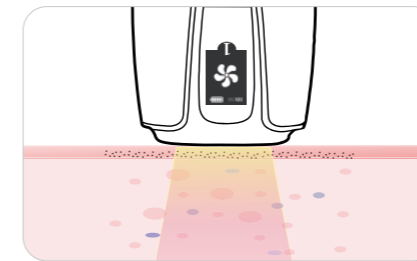

The (IPL) intense pulsed light stimulates regeneration of collagen and elastic fiber thereby improves the color and texture of the skin without a surgery

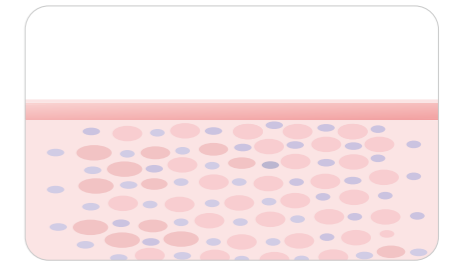


The IPL can can treat some of the visible damage caused by the sun exposure called photoaging, which you may notice mostly on your face,neck,hands,or chest. After using IPL your skin tone gets more even, reduction of wrinkles and skin rejuvenation.


Contraindications

1. General Conditions

- ❖ Never use the device if you have dark skin. Check the skin tone table below to see if your skin allows usage of the device.

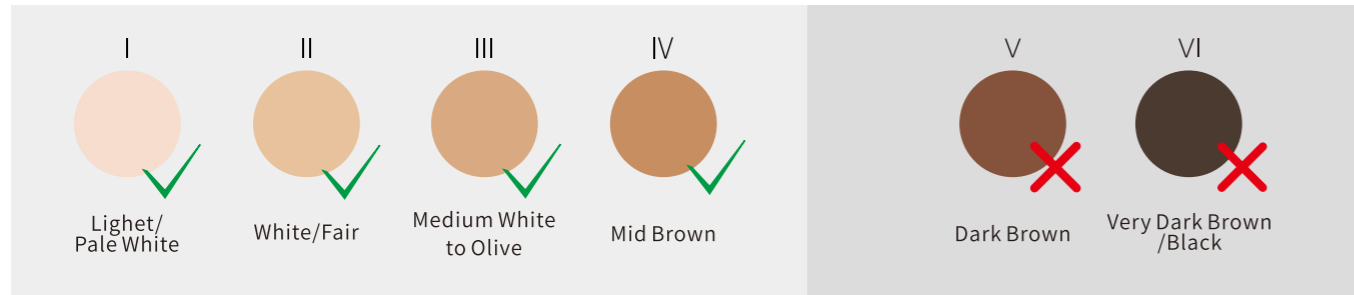

- ❖ The device is not effective if your hair color is light brown, light gold, red or white gray. Check the hair color table below to see if your hair color allows usage of the device.

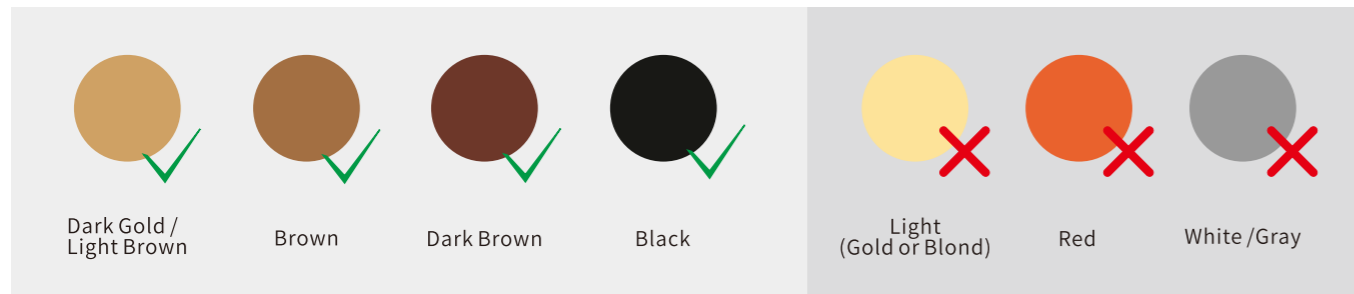

- ❖ Never use the device if you are pregnant or breastfeeding.
- ❖ Never use the device if you have active implants.

Contraindications

2. Medications / History

- ❖ **Never use the device if you are taking any of the medications/treatment listed below:**

I) If your skin is currently being treated with or has recently been treated in the past week with Alpha-Hydroxy Acids (AHAs), Beta-Hydroxy Acids (BHAs), topical isotretinoin and azelaic acid.

II) If you have taken any form of isotretinoin Accutane or Roaccutane in the last 6 months. This treatment can make skin more susceptible to tears, wounds and irritations.

III) If you are taking photosensitizing agents or medications, check the package insert of your medicine and never use the device if it is stated that it can cause photo-allergic reactions, photo-toxic reactions or if you have to avoid sun when taking this medicine.

IV) If you take anticoagulation medicine, including heavy use of aspirin in a manner which does not allow for a minimum 1-week washout period prior to each treatment. White / Fair Medium White to Olive Mid Brown Dark Brown Very Dark Brown / Black Dark Gold / Light Brown Brown Dark Brown Black Light (Gold or Blond) Red White / Gray

V) If you have received radiation therapy or chemotherapy within the past 3 months.

VI) If you are on painkillers which reduces the skin sensitivity to heat.

VII) If you take immunosuppressive medications.

VIII) If you have or had surgery in the areas to be treated for IPL in the last 3 weeks.

3. Disorders

- ❖ **Never use the device;**

I) If you have diabetes or other systemic or metabolic diseases.

II) If you have congestive heart disease,

III) If you have a disease related to photosensitivity.

IV) If you have a history of collagen disorder.

V) If your skin is sensitive to light and easily develops a rash or an allergic reaction.

VI) If you have a skin disease such as active skin cancer or you have history of skin cancer or any other localized cancer in the areas to be treated by IPL.

VII) If you have a history of vascular disorder.

VIII) If you have any bleeding disorder.

IX) If you have a history of immunosuppressive disease

Contraindications

4. Skin Condition

❖ Never use the device;

- I) If you have infections,eczema,burns,inflammation of hairfollicles,open lacerations,herpes simplex (cold sores), wounds or lesions and hematomas in the area to be treated by the IPL.
- II) On irritated (red or cut),sunburned,recently tanned or fake-tanned skin.
- III) On these areas (without consulting your doctor first) such as;moles,freckles,large veins,darker pigmented areas,scars and skin anomalies. This can result in a burn and a change in skin color,which makes it potentially harder to identify skin-related diseases.
- IV) On areas such as; warts,tattoos or permanent makeup.

5. Areas

❖ Never use the device;

- I) Around the eyes or near the eye brows.
- II) On the lips,nipples,areola,labia minora,vagina,anus and inside of the nostrils and ears.
- III) Men must not use the device on the face and neck including all beardgrowing areas,nor on the whole genital area.
- IV) On areas where you use long lasting deodorants.
- V) Over or near anything artificial like silicone implants,subcutaneous injection ports (for instance an insulin dispenser) or piercing.

⚠ **If you are not sure whether you can use the device, please consult your doctor.**

Important Notes

1.Warning

- ❖ Keep the device and the adapter dry.
- ❖ If the device is broken, do not touch any inner part to avoid electric shock.
- ❖ Water and electricity are a dangerous combination, Do not use this device in a wet surroundings.
- ❖ This device is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, unless they have been given supervision or instruction concerning use of the device by a person responsible for their safety.
- ❖ The device is not intended for children under the age of 15 years.
- ❖ Always check the device before you use it. Do not use the device or adapter if it is damaged .
- ❖ Do not use the device, if the flash(light) exit window is broken.
- ❖ The flash(light) exit window can become very hot after use. Do not touch the inner part of the light exit window, before it cools down.

2.Device Care

- ❖ Make sure that nothing obstructs the airflow through the ventilation slots of the device.
- ❖ Store the device in a dust free and dry place.
- ❖ To prevent damage, do not expose the device to direct sunlight or UV light for several hours.

3.Caution

- ❖ For hygienic reasons, the device should only be used by one person.
- ❖ Use the device only at settings suitable for your skin type. Using at higher settings than those recommended can increase the risk of skin reactions and side effects.
- ❖ Do not look directly at the flash window while using the device. It is not necessary to wear goggles however as a precaution goggles are provided along with the device for your safety.
- ❖ Use the device in a well-lit room, so that the light is less glaring to your eyes.
- ❖ Never leave the device unattended when it is switched ON. Always switch OFF the device after use.
- ❖ Tanning with natural or artificial sunlight might influence the sensitivity and color of your skin.
- ❖ Before you use the device, you should clean your skin and make sure it is hair-free,entirely dry and free from oily substances.
- ❖ Do not treat the same skin area more than once during a session.
- ❖ The device should never be painful. If you experience discomfort, reduce the light intensity level.

Course Of Treatment

1st Treatment - (Week - 0)

- ❖ After the first treatment, it can take 1 to 2 weeks for the hairs to fall out.
- ❖ During the initial treatments, you might still see some hairs growing. These are likely to be the hairs that were in their growing phase during the first treatment.

2-3 Treatment - (Week - 2 & Week 4)

- ❖ After 2-3 treatments, you should see a noticeable reduction in hair growth. However, to effectively treat all hairs, it is important to keep on treating according to the recommended treatment schedule.

4-5 Treatment - (Week - 6 & Week 8)

- ❖ After 4-5 treatments, you should see a significant reduction of hair growth in the areas that you treated with the device.
- ❖ A reduction of hair density should also be visible. Keep on treating with frequent touch-ups (every 4 weeks) to maintain the result.


Once every two weeks

Once a month


Use when needed




Zero week

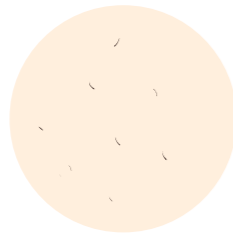

Two weeks

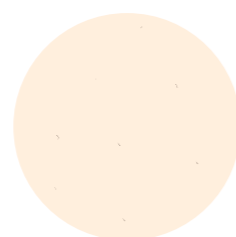

Four weeks


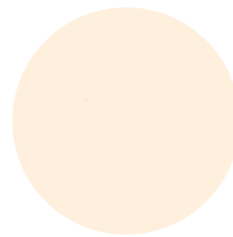


Six weeks

Before Use

1. Get your skin ready

- ❖ Before you use the IPL device, remove all visible hair in the area you wish to treat.
- ❖ Ensure there are no hairs remaining above the surface of the skin. This allows the light to be easily absorbed by the melanin located in the hair roots and ensures effective treatment.
- ❖ You should only shave before using the IPL.
- ❖ Do not wax, epilate or depilatory creams, as anything else pulls hair from the roots (follicle) and hinders the IPL process, which you don't want.
- ❖ If you choose to wax, please wait 6-8 days before using the device to let your skin rest.
- ❖ We recommend that you take a shower before the treatment to ensure the skin is free from oily substances, dirt and then pat your skin dry.

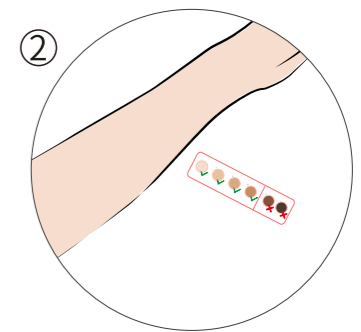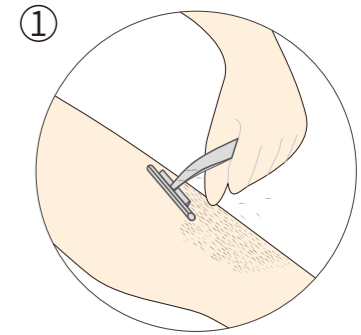


2. Skin test

- ❖ Check the skin tone table to see if your skin type allows the usage of the device.



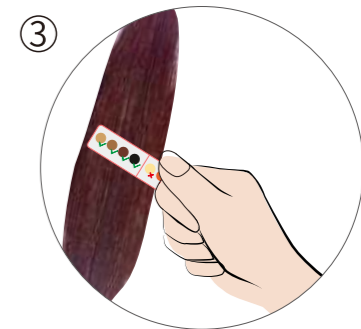

3. Hair test



- ❖ Check the hair color table to see if your hair color allows the usage of the device.



While Use


1. Selecting the right intensity

- ❖ The device provides 8 different light intensities and advises you on the right setting you should use based on your skin tone. You will always be able to select the light intensity settings as you find comfortable.
- ❖ Press  for 2 seconds to turn the device ON
- ❖ Short press  to adjust the light intensity level, Each time you press  the light intensity level will gradually Increase from 1st gear to 8th gear. After 8th gear, press  the light intensity level will gradually decrease until 1st gear.

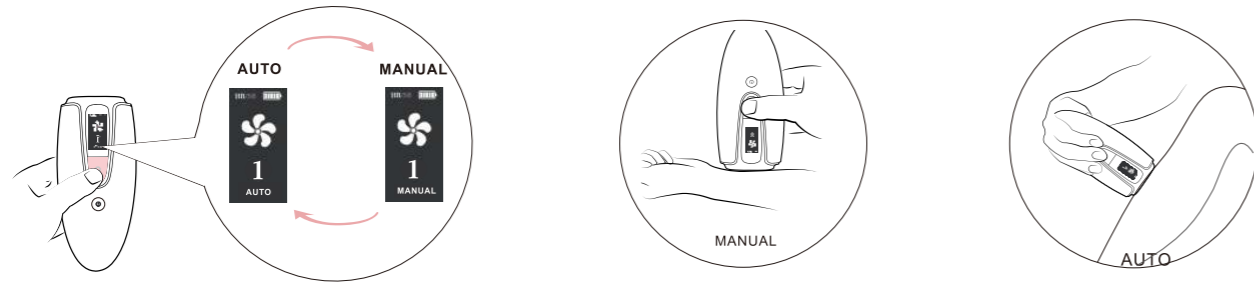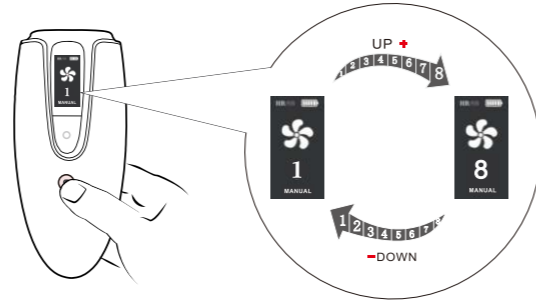

- ❖ Press  and  hold for 2 seconds to switch mode between "Auto" and "Manual"

Note: There is only one intensity level for skin rejuvenation.

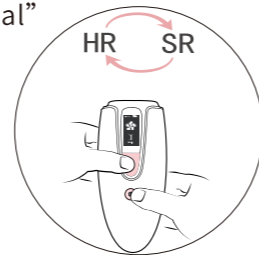


While Use

For hair removal, select the desired light intensity settings. Please see the table below to set the light intensity settings.


Skin Type & Tone	Light intensity Setting
I White: always sunburn, never tan.	6-8
II Beige: easily sunburn, tan minimally.	5-7
III Light brown: sunburn easily, tan slowly to light brown.	4-6
IV Middle brown: rarely sunburn, tan easily.	2-4
V Dark brown: rarely sunburn, tan very easily.	You cannot use the device
VI Brownish black or darker: rarely or never sunburn, very dark tanning.	You cannot use the device

- ❖ The IPL should never be painful. If you experience discomfort, reduce the light intensity settings.
- ❖ After moving to another body area or recent tanning, perform a skin test to determine the right light intensity settings.


2. Handling the device

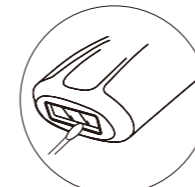

1. Clean the light exit window before use

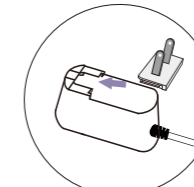


2. Insert the plug

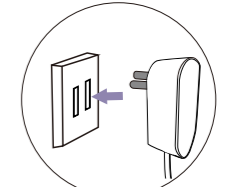


3. Plug the power supply in the wall socket



4. Switch on the device and select the right light intensity for your skin tone

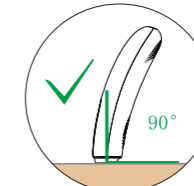

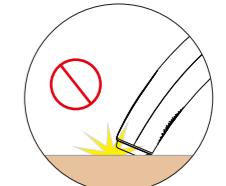

5. Place the device at a 90° angle on the skin so that the skin contact surface of the device is in contact with your skin



6. leaving no gaps between skin and device in case of blinking


While Use

- Press the device firmly onto your skin to ensure proper skin contact. Press  to release a flash. You will feel a warm sensation because of the flash and its normal.
- Place the device (flash/light) window in full contact with the skin to make sure the product works correctly. If not in proper contact with the skin, the device will not work.
- Place the device on the next area to be treated. After each flash, it takes up to 1-3 seconds until the device is ready to flash again.
- For larger body areas use AUTO mode. Hold down the button and treat along your skin. For smaller areas use MANUAL mode.
- When you have finished the treatment, press and hold  for 2 seconds to switch OFF the device.
- Remove the adapter from the electric socket.
- It is important to wipe clean the (flash/light) window after every use.

3. Treatment Schedule

1>. Hair removal schedule

- ① In the first phase of your treatment, treat once a week for 4 - 8 weeks. (Depending on your skin result)
- ② After this period use your handset once a month, for 3 months or until satisfied.
- ③ To maintain hairless smooth skin, use your handset once every 2 to 3 months, or as needed.

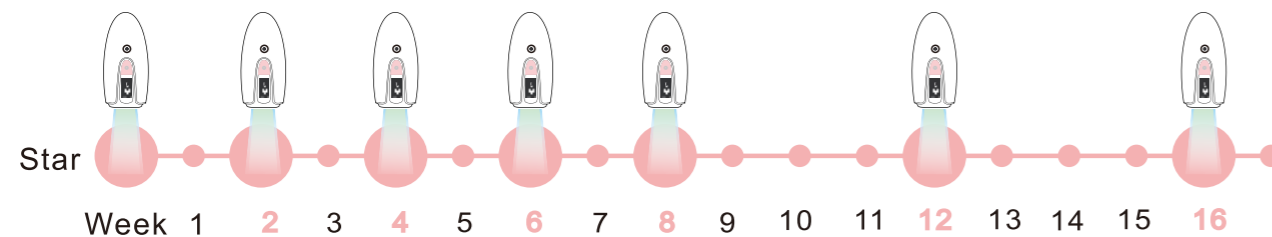

2>. Skin rejuvenation schedule

Please note that for skin rejuvenation, we suggest to use every 3 days.

After Use

1. After Care

- ❖ After the IPL treatment, you can safely apply lotion, creams or moisturizer on the treated area.
- ❖ If you experience skin irritation or skin redness after the treatment, wait until it disappears before applying any product (lotion or cream) on your skin.
- ❖ If you experience skin irritation or skin redness after applying any product (lotion or cream) on your skin, wash it off with water.

2. Common reactions

- ❖ If your skin shows slight redness and/or may prickle, tingle or feel warm. (The reaction is absolutely harmless and will disappear soon).
- ❖ Dry skin and itching may occur due to shaving or a combination of shaving and IPL treatment. (You can cool the area with an ice or a wet face cloth. If dryness persists, you can apply a non-scented moisturizer on the treated area.)

3. Rare Side effects

- ❖ Burns, excessive redness (e.g., around hair follicle).
- ❖ Skin discoloration: this occurs very rarely.
- ❖ Epidermal heating

Device Cleaning And Storage

- ❖ After usage - Switch OFF the device , unplug it and let it cool down.
- ❖ Use a soft, moist cloth to clean the IPL parts such as; (flash/light) window, the outside surface of the attachments and the reflectors on the side.
- ❖ Store the device in a dust free and dry place at a temperature between -25°C and 75°C

Device Specifications

Product	IPL hair removal device	Spectrum	>550nm
Model	DH5-IPL	Pulse duration	0.9+/-0.2ms
Rated input	AC100-240V/50-60Hz	Max light energy	7.92-15.31J
Rated output	DC 12V, 2.5A	Operating Temperature	5-35C
Executive standard	GH4706.1 GB4706.15	Operating Humidity	15%-75%
Protection against electric shock	II	Weight	209g
Protections rating	IP 30(EN 60529)	Size	H179*W70*D35mm



Frequently Asked Questions

Q. The IPL device does not flash,when I place the device on the skin?

A: Please make sure the device is at a 90 degree angle on the skin,so the device sensor can detect the skin and thereby flash.

Q. Why the device produces strange smell ?

A: Either the device light/flash (window) was dirty or you did not shave the treatment area properly. So the burned hair produced the strange smell.

Q. Why my skin reaction after treatment last longer than usual ?

A: You may have used an IPL intensity level settings that is too high for your skintype.

Q. I don't see any significant results after the treatment?

A: You may have used an IPL intensity level settings that is too low for your skin type OR you respond more slowly towards your IPL treatment.

Q. Why does the device/adapter becomes warm during use?

A: This is normal. The device produces large amount of energy when it emits light,which causes the device to warm-up

Q. Why should we shave before using this IPL device?

A: The goal is to reach the IPL light energy to the hair follicle as quickly as possible. So shaving before is essential to achieve the best results.

Frequently Asked Questions

Q. Can I use the IPL device before or after tanning?

A: If you have recently been tanning, we recommend waiting 2 weeks before using your handset. You should wait 48 hours to tan after using your handset. Artificial tan should not be on your skin while using your handset.

Q. Is this suitable for men?

A: Definitely! Our handsets are suitable for both men and women.

Q. How often will I have to treat to maintain the hair removal results?

A: After completing your start-up phase, you can continue to treat according to your needs. Retreat on a per need basis.

Q. Do I need to stay out of the sun after using IPL device?

A: After treatment, we recommend using sunscreen SPF 15, whenever treated areas may be exposed to the sun.

Q. My skin turned red and feel pain after using the IPL device. Is this a problem?

A: Treated skin may temporarily turn red or feel painful depending on the skin type or on the method the IPL device is used by the person.If the above symptoms do not disappear in 2 days ,consult a dermatologist.

Device Warranty Card

Owner Warranty Certificate

(User save, show it to the maintenance point of repair)

Product	IPLHair Removal Device	No.	
Warranty	1 Year warranty start from the date of purchase		
Name		Tel?	
Social	Wechat / QQ	facebook	
account Address ?			
Date of purchase		Warranty Deadline	
Purchased from	Shop?		
	Address?		
	Tel?		
Purchase sheet	InvoiceNo.	Order No.	
Problem Description			
This product is only used for household usage			