



USER MANUAL

Titan Gold



Table of contents

1.	Introduction -----	P.3
2.	Caution -----	P.3
3.	Warning -----	P.4
4.	Possible Side Effect -----	P.6
5.	Precaution of Safe Use -----	P.9
6.	Presentation of Your IPL Hair Removal System	
6.1	Working Principle -----	P.13
6.2	The Impact of Growth Cycle of the Hair -----	P.14
6.3	Select the Intensity Level with Caution -----	P.15
6.4	Session Planning -----	P.17
6.5	Expected Epilation Result -----	P.18
6.6	Using the IPL Hair Removal System Safely -----	P.18
7.	Use -----	P.19
8.	Maintenance of Your IPL Hair Removal System -----	P.26
9.	Troubleshooting -----	P.27
10.	Customer Service -----	P.28
11.	FAQ -----	P.29
12.	Specification -----	P.33

1. Introduction

The IPL Hair Removal is a powerful electrical appliance and so for reasons of safety it should be used with great care.

Before using the Hair Removal for the first time, please carefully read this user manual completely and pay attention in particular to the safety instructions and the precautions of use.

2. Caution

2.1 This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance.

2.2 Cleaning and user maintenance shall not be made by children without supervision.

2.3 The appliance is only to be used with the power supply unit provided with the appliance.

2.4 The appliance is intended for households and indoor use only.

2.5 Appliances not suitable for bath or shower use.



2.6 A warning that intense visible and invisible optical radiation emitted from the device may cause eye injury (IEC 60335-2-113)

2.7 It is strictly forbidden to expose this equipment in the eye area;

2.8 It is strictly prohibited to play with children or infants as toys. When unused, they are not available in children.

2.9 Warning against the use of the device may cause eye (or even blindness) or skin injury or discomfort if it is not operated according to the instructions. Please follow the instruction manual when using.

2.10 Warning prolonged use or repeated use of the same skin may cause skin damage

2.11 Incorrect use of the device may cause eye damage, and it is strictly prohibited to use the device around the eye or around the eye.

2.12 If you are taking a photosensitive medicine or herbal medicine, consult your doctor before using the device to avoid injury such as skin damage

2.13 It is prohibited to use the device in areas where the skin is damaged, such as burns, injuries, blister or edema, scarring, wound, skin disease, sunburn, infection, etc

2.14 Warning the device is not suitable for therapeutic purposes. Consult your doctor for assistance, such as hemorrhoids, skin rashes, skin itching, skin fungus, skin redness and skin irritation, etc.

2.15 Warning that the device does not apply to children or infants;

2.16 It is not possible to use the device if it is not visible to the part of the equipment that has been used to confirm the luminescence of the equipment or the launch box.

2.17 Warning when using the equipment, cannot use lotion or cosmetics, banning the use of flammable liquid, such as alcohol preparation (including perfume, deodorant, or any other alcoholic) or. Acetone, clean your skin, and skin surface when using fur to shave in advance.

2.18 Note that the equipment can be cleaned with alcohol or dry strips and placed in a dry, dry 2.19 Warning to consult your doctor for advice when using the device for a cycle.

2.20 White or gray hair is not suitable for this device.

2.21 Warning that minors should be under parental supervision to clean and maintain the equipment.

2.22 Warning over the age of 8 minors and not using experience, have a physical or mental handicap, light sensitive person must understand the premise of product hazards, to operate the equipment safely under supervision and guidance.

3. Warning

Before using your IPL Hair Removal appliance, please read thoroughly and follow precisely all safety warnings .

Note: The IPL Hair Removal appliance is not as effective on white or light – colored hair. We recommend that you read this user manual before every use of your IPL Hair Removal. If you are unsure about using the IPL Hair Removal , please consult your dermatologist or your GP.

3.1 Do not use the IPL Hair Removal on naturally dark skin!

3.2 Do not use the Hair Removal appliance on tanned skin or after recent exposure to the sun! Avoid exposing yourself to sunlight for 4 weeks before treatment with the IPL Hair Removal, and for 2 weeks following treatment.

3.3 Do not use the IPL Hair Removal on white, grey or blond hair.

3.4 Do use the IPL removal system on the scalp, eyebrows, chin, near the eyes, ears, nose, nipples,

genitals, anus, as well as in areas with various marks, freckles, moles, warts, bruises, rashes, tattoos or permanent make-up, etc.

3.5 Do not use if you are pregnant or breastfeeding

3.6 You have had lengthy exposure to the sun or UV equipment over the previous 28 days

3.7 Do not use it if you suffer from eczema, psoriasis, lesions, open wounds or evolutive infection in the area to be treated. Wait until the affected area has healed before using the IPL Hair Removal Appliance.

3.8 Do not use if You show antecedents to the genesis of keloids, you have a known sensitivity to light (photosensitivity) or one caused by medication, including non-steroidal anti-inflammatories (e.g.: aspirin, ibuprofen, paracetamol), tetracycline, phenothiazine, thiazide, diuretics, sulphonylureas, sulfamides, DTIC, fluorouracil, vinblastine, griseofulvin, Alpha Hydroxy Acids (AHAs), Beta-Hydroxy Acids (BHAs), Retin-A R, Accutane R or topical retinoids.

3.9 Do not use it if you have recently been treated or being treated with Alpha-Hydroxy Acids (AHA), Beta-Hydroxy Acids (BHA), Retin-A R, topical retinoids or azelaic acid.

3.10 Do not use it if you suffer from an abnormal skin condition caused by diabetes, for example, or another systemic or metabolic disease.

3.11 Do not use it if you have been undergoing treatment with Accutane R (isotretinoin) for the previous 6 months.

3.12 Do not use it if you have been following a steroid regime for the previous 3 months.

3.13 Do not use it if you show signs of pre-herpes growth in the area to be treated, unless you have consulted a physician and been given preventive treatment before using the IPL Hair Removal .

3.14 Do not use it if you suffer from epilepsy.

3.15 Do not use it if you have an active implant, like a pacemaker, an implant for urinary incontinence, an insulin pump, etc.

3.16 Do not use it if you have a pathology linked to photosensitivity, like porphyria, polymorphic eruption to light, solar urticaria, lupus, etc.

3.17 Do not use it if you have an antecedent to skin cancer or areas subject to dermal malignancies.

3.18 Do not use it if you have been undergoing radio or chemotherapy for the last 3 months.

3.19 Do not use it if you suffer from any other condition your physician deems a risk for the treatment.

3.20 Wait until the affected area has healed before using the IPL Hair Removal appliance.

If you have any doubts about the risks to you of the use of the IPL Hair Removal appliance, consult your GP or your dermatologist.

4. Possible side effects

When your IPL Hair Removal is used according to these instructions, side effects and complications are rare. Nonetheless, all cosmetic treatments, even those designed for use at home, involve some risks. Therefore, it is important to understand and accept these risks and possible complications when you use a pulsed light depilation system designed for use at home.

If you superimpose pulses or administer multiple pulses on the same patch:

Side Effects	Appearance of the Side Effect
Slight discomfort/annoying sensation in the skin	light
Redness	light
Hypersensitivity of the skin	light
Burns or sores	Rare
Scarring	Rare
Change in pigmentation	Rare
Excessive redness of the skin and swelling	Rare
Infection	Negligible
Bruising	Negligible

4.1 Slight annoying sensation in the skin.

Even if the IPL Hair Removal depilation is generally tolerated well, most people who use it feel a slight annoying sensation during use, often described as a prickling sensation on the area being treated. Generally, this sensation lasts throughout the application. Even for a few minutes afterwards.

Any discomfort beyond this annoying sensation is abnormal and means that either you must stop the treatment , because you cannot tolerate it, or you must lower the selected power level.

4.2 Redness

Your skin could become red just after use of the IPL Hair Removal or in the 24 hours following treatment. The redness generally disappears within 24 hours. Consult your physician if it persists longer than 2 to 3 days.

4.3 Hypersensitivity of the skin

The skin is more sensitive in the area that has been treated and could dry out and scale.

4.4 Burns or sores.

Very rarely the skin may be exposed to injuries or burns following the application. Several weeks may be necessary for it to heal completely and, in exceptional cases, a scar may remain.

4.5 Scarring

Even if these cases are very rare, permanent scars may appear. Generally the scars take the form of a whitish superficial lesion (hypotrophic). Nonetheless, the scar may be large and red (hypertrophic) or large and extended beyond the injured area (keloid). More advanced cosmetic treatment may be required to improve the appearance.

4.6 Change in pigmentation

The IPL Hair Removal appliance targets the hair roots, in particular the pigmented cells of the bulb. Nonetheless, there is a risk of temporary hyperpigmentation (increase in the brown pigmentation or discoloration) or Hypopigmentation (Whitening) of the surrounding skin. The risk of change in pigmentation is higher among people with darker skin.

Generally, the change in pigmentation and the discolouration are temporary. Cases of permanent hyperpigmentation or hypopigmentation are rare.

4.7 Excessive redness and swelling

In some rare cases, treated skin may become very red and swollen. This reaction is more common in sensitive areas of the body. The redness and the swelling should disappear within 2 to 7 days during which you should apply ice to the affected area regularly. Gentle cleaning does not pose a problem, but you should avoid any exposure to the sun.

4.8 Infection.

A skin infection is extremely rare, but may occur after a burn or a lesion of the skin caused by the IPL Hair Removal

4.9 Bruising

Very rarely, using the IPL Hair Removal may cause some blue or purple bruising which could last 5 to 10 days. As the blue disappears, the skin may take on a rusty colour (hyperpigmentation) that could be permanent.

Please consult your doctor if you experience persistent discomfort or irritation during or after the IPL Hair Removal treatment.

5. Precautions for safe use

Follow the instructions in this manual very carefully when selecting the power intensity during each session with the IPL Hair Removal

Selecting the inappropriate power intensity may cause side effects.

For detailed instructions on selecting the power, see the item "select the power intensity level with "CAUTION" in the section "Using the IPL Hair Removal "

5.1 Avoid side effects!

Do not treat the same patch of skin more than once per session!

Avoid using a pulse on the same patch!

If burns or blisters appear, STOP IMMEDIATELY!

5.2 Avoid complications after using the IPL Hair Removal !

Avoid exposure to sunlight. "Exposure to sunlight" means continuous exposure to direct sunlight without protection for about fifteen minutes or continuous exposure to indirect sunlight without protection for an hour;

Avoid exposure to the sun for 4 weeks prior to a treatment with the IPL Hair Removal and again for 2 weeks following the treatment.

Appropriate protection, like clothing that offers sufficient coverage or a sunscreen (with a protection factor higher than 50). Should be used regularly 3 weeks before and after the treatment.

Caution: It is important that you follow these instructions carefully if there is exposure to UV on a sun bed.

- If you have applied a self- tanning cream to your skin, wait until the cream's tanning effect has disappeared completely before using your appliance.
- Always shave the area to be treated and make sure you have clean skin before using the IPL Hair Removal appliance.
- Make sure birthmarks and tattoos are covered during the IPL Hair Removal treatment and also cover brown or black marks, as well as freckles, nevus, warts or growths before using the IPL Hair Removal
- Never look directly into the light emitted by the flash lamp. Use the safety glasses included whenever you use the IPL Hair Removal .
- Never use the IPL Hair Removal on the tip of the breasts or the genital organs (male or female).
- Never use the IPL Hair Removal on any part of the body where you might later want hair to grow back.
- Never use the IPL Hair Removal for anything but depilation.
- Never point the application surface of the IPL Hair Removal to emit a light pulse into open space.
- Always ensure you point the treatment surface at the skin and that it is in full contact with your skin before releasing a pulse.
- Caution: surfaces are hot. Do not touch the lens on the applicator surface with your finger. Also ensure you keep the appliance's surface clean. See point 7.1 on cleaning and maintenance of your IPL Hair Removal .
- Never use flammable liquids, like alcohol (including perfume, deodorant or any other preparation containing alcohol) or acetone, to clean you skin before using the IPL Hair Removal Using your IPL Hair Removal appliance may cause a temporary change in pigmentation (see below Possible Side Effects).
- Keep the appliance out of the reach of children. Do not use the IPL Hair Removal on children nor allow children to use it.
- This appliance is not intended for use by people(adults and children) whose physical, sensory or mental capacities are reduced, or by persons who lack experience or knowledge, unless they are supervised or have been given instructions beforehand on the use of the appliance by a person responsible for their safety . Children should be supervised to ensure that they do not use the appliance like a toy.



- **CAUTION!:** The plastic bags used to wrap this appliance or the packaging may be dangerous. Keep these bags out of the reach of babies and children Do not use them in cribs, children's beds, pushchairs or playpens. The thin film will stick to the nose and mouth and prevent breathing. A bag is not a toy.



- If using the appliance in a bathroom, make sure you unplug it after use. In fact, having it close to a water source may be dangerous even if the appliance is turned off . To ensure additional protection, we recommend installing a residual current device (RCD) with a

nominal trip current of no higher than 30 mA in the electrical circuit that supplies power to the bathroom. Ask your electrical fitter for advice.

- ELECTRICAL AND ELECTRONIC EQUIPMENT AT THE END OF LIFE.
- In everybody's interest and to participate actively in protecting the environment:



- Do not discard these products with your household waste.
- Use return and collection systems available in your country.
- Some materials can in this way be recycled or recovered.
- Do not use the unit if it has fallen or if it appears to be damaged.
- Do not use any accessories other than those recommended by factory

As with all electrical appliances, some precautions are necessary to ensure your safety when you use the IPL Hair Removal .



- Keep the IPL Hair Removal away from all sources of water!
- The IPL Hair Removal is an electrical appliance and should never come in contact with water.
- Do not store the IPL Hair Removal in a place where it could drop or fall into a bath or washbasin or any other container with liquid in it. Do not put it in water or any other liquid. This could result in severe electrical shock.
- Do not use your IPL Hair Removal when having a bath.
- Never use the IPL Hair Removal appliance if it is damp or wet.
- Do not catch hold of the IPL Hair Removal if it has fallen in water.
- Unplug the IPL Hair Removal appliance immediately if it has fallen in water.



- Do not attempt to open or to repair your IPL Hair Removal appliance yourself . Opening the appliance will expose you to dangerous electrical components and to the energy of the pulsed light, which could cause serious physical injury and /or irreversible eye damage.

Only repair centers certified by the factory are authorised to repair any damage.

- Please contact the customer service of your IPL Hair Removal if our appliance is broken, damaged or needs repair.
- Only use the IPL Hair Removal for what it was designed to do and following the instructions of this manual.
- Your appliance should never be left unsupervised while plugged in.
- Never leave your IPL Hair Removal unsupervised while turned on.
- Do not plug your IPL Hair Removal in if the cord or the plug is damaged. Keep the cord away from all heat sources.
- Stop using your appliance if you notice a defect or if it appears defective or emits an abnormal noise or odour .In particular, stop using the appliance immediately if the flash lamp is cracked or is stained.
- Stop using the appliance when the READY indicator light turns red and disconnect the IPL Hair Removal appliance adaptor from the mains outlet.
- Always disconnect the IPL Hair Removal appliance adapter from the mains outlet immediately after use.
- Unplug your IPL Hair Removal appliance before you clean it.
- Use only attachments and parts with our IPL Hair Removal appliance that have been recommended by Factory.
- Only use the adaptor provided to operate the unit.

6. Presentation of your IPL Hair Removal System

6.1 Working Principle

For many years, the factory has placed itself at the leading edge of technology to better respond to the increasingly demanding and pointed desires of women in depilation. Thanks to the IPL(Intense Pulsed Light) epilation technique, the IPL Hair Removal will eliminate unwanted hair gently and painlessly leaving smooth and silky skin .

The IPL Hair Removal is a personal long-lasting epilation device. The IPL Hair Removal procedure is

well known and has been mastered. It has been subjected to over 15 years of clinical testing around the world and has proven its effectiveness and its safety in offering gradual but long-lasting hair removal.

IPL Hair Removal based on the theory of selective photothermolysis which uses optical energy to stop the growth of the hair. To obtain such thermal effects, the root of the hair must selectively absorb the light energy and transform it into heat. This selection is produced when most of the optical energy emitted is absorbed by the pigment of the bulb, while the surrounding skin and tissue remain cold. This thermal effect painlessly destroys the hair and roots and prevents regrowth. This scientifically proven technology gives quick and lasting results.

The amount of pigment in a person's skin determines not only the colour, but also the level of risk to which the user is exposed when using a form of depilation that uses light, regardless of which one. The darker the skin, the higher the risk will be.




6.2 What impact does the growth cycle of the hair have on photo-epilation

Melanin is the pigment found in the hair, which allows the absorption of the light, to transform it into heat in order to stop the hair growth. Also, the more melanin in the hair (i.e., the darker the hair), the more light it can absorb and the more efficient photo-epilation becomes.

Every hair on our body has three growth phases: the anagen, catagen and telogen phases. These three phases play an important role in the process of IPL-epilation devices. The anagen phase is the growth phase of the hair, while the catagen and telogen phases are both resting phases. The time the hair takes to go from one part of the body to the next. We can, however, establish an average of between 18 and 24 months.

Most hair follicles are found at any one time in one of the two resting phases. Consequently, the IPL Hair Removal will have no effect when they are in these resting phases. Nonetheless, it will affect the hair that is in anagen phase. It is therefore important to understand the duration of a full follicular cycle to obtain lasting results with your IPL Hair Removal.

The three growth phases: the anagen, catagen and telogen phases

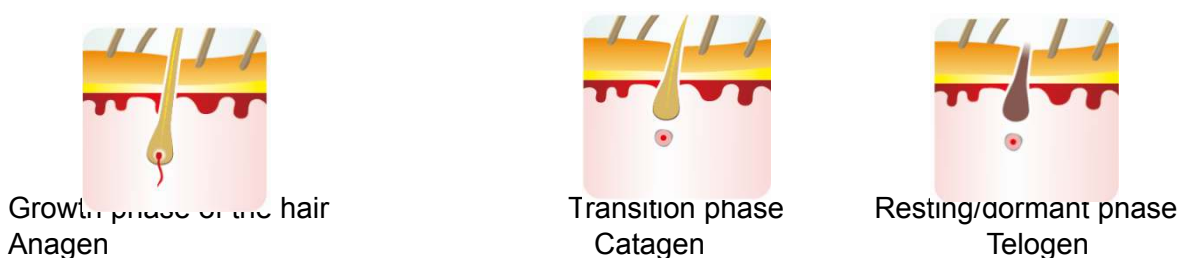

6.3 Select the intensity level with CAUTION

The power intensity level designates the power of pulsed light administered to your skin during treatment ranging from the lowest (1) to the highest (5) levels.

The power indicator lights indicate the selected level of power. If you increase the power the results will be better, but there will also be an increased risk of side effects (cf. the section "possible side Effects")

CAUTION: 48 hours correctly to the test, i.e. if you do not notice any abnormal reactions or discomfort in the 48 hours following the test, you can select the power level of your appliance for your first session using the following chart:

- For the first session, people with phototype I or II skin should set the appliance to level 3 at most. Gradually, over the next few sessions, if you do not experience any specific discomfort or unpleasantness, the power level can be increased to level 4 and then to level 5.
- For the first session, people with phototype III or IV skin should set the appliance to level 2 at

most. Gradually, over the next few sessions, if you do not experience any specific discomfort or unpleasantness, the power level can be increased to level 3 then to level 4.

- For the first session, people with phototype V skin should set the appliance to level 1 at most. Gradually over the next few sessions, if you do not experience any specific discomfort or unpleasantness, the power level can be increased to level 2 then to level 3.
- The IPL Hair Removal appliance is not suitable for use by people with phototype VI skin.



Read the instructions fully before using the appliance and refer to the table below representing various phototypes to select the safe power level every time you carry out a treatment.

The intensity selection guide must be used along with the table below for a safe use of the IPL Hair Removal appliance.

PHOTOTYPES (according to the Fitzpatrick scale):

- I: Very white, milky skin- Never tans and always burns
- II: Light – coloured skin – Tans lightly and usually burns
- III: Beige skin – Tans moderately and sometimes burns
- IV: Medium brown skin – Tans easily and seldom burns
- V: Dark brown skin – Tans very easily and rarely burns
- VI: Black skin

SKIN CHART

Skin Type	Skin Colour	Skin Tone	Tanning History	Maximum Power Level
I	Very fair, 'transparent'		Always burns, never tans	Subject to patch test, all 5 power levels can be used
II	Fair		Always burns	
III	Fair to light olive		Burns mildly	Subject to patch test, power levels 1, 2, 3 and 4 can be used
IV	Olive/Tanned		Rarely burns, tans with ease	
V	Dark brown		Very rarely burns	Subject to patch test, power levels 1, 2 and 3 can be used
VI	Black		Never burns	NOT suitable for IPL treatment

PHOTOTYPES	I	II	III	IV	V	VI
SESSIONS Test (carried out 48 hrs before the first session)	Level 1	Level 1	Level	Level 1	Level 1	X
First session	Level 3	Level 3	Level 2	Level 2	Level 1	X
Following sessions	Level 4-5	Level 4-5	Level 3-4	Level 3-4	Level 2-3	X

Note: X--The appliance is not suitable for use on phototype VI skin

It is possible that the power intensity level required will vary with the difference in a person's skin colour to the area to be treated. Skin on the arms, generally exposed to the sun more frequently, is for example darker than the underarm.

But be careful by increasing the power level for light areas on your body. You will have better results but you also increase the risk of discomfort or even side effects.

6.4 Session planning

A full follicular cycle takes on average from 18 to 24 months during which many sessions with your IPL Hair Removal will be necessary to obtain long-lasting results. The efficiency of the appliance will vary from one person to the next and will depend on the area being treated, the colour of the hair and the way in which the appliance is used.

Below is an example of a depilation schedule over a full growth cycle.

- The first 4 sessions with the IPL Hair Removal should be 2 weeks apart.
- Next sessions should be 4 weeks apart until you achieve the desired results.

1st month						
1	2	3	4	5	6	7
8	9	10	11	12	13	14
15	16	17	18	19	20	21
22	23	24	25	26	27	28
29	30					
2nd month						
1	2	3	4	5	6	7
8	9	10	11	12	13	14
15	16	17	18	19	20	21
22	23	24	25	26	27	28
29	30					

3rd month						
1	2	3	4	5	6	7
8	9	10	11	12	13	14
15	16	17	18	19	20	21
22	23	24	25	26	27	28
29	30					
4th month						
1	2	3	4	5	6	7
8	9	10	11	12	13	14
15	16	17	18	19	20	21
22	23	24	25	26	27	28
29	30					

5th month						
1	2	3	4	5	6	7

8	9	10	11	12	13	14
15	16	17	18	19	20	21
22	23	24	25	26	27	28
29	30					

After that you need only use the IPL Hair Removal when hair regrows until epilation is lasting. Leave at least 4 weeks between sessions.

You can plan your personal sessions to suit your rhythm and still achieve satisfactory results. However, we recommend you follow the schedule suggested here because the results will be excellent.

Note: Do not treat the same patch more than once in two weeks as this will not improve results and it will risk causing undesirable effects.

You can plan maintenance depilation sessions with your IPL Hair Removal once you have achieved your objectives. Because of hormonal factors and physiological changes, dormant follicles can become active. Maintenance sessions may thus sometimes be necessary.

6.5 Expected epilation results

The first results are visible after 4 sessions and achieve a very significant decrease in the amount of hair after 9 to 10 sessions. The result will vary from one person to the next and according to the kind of hair and biological factors that determine the growth cycle of the hair. Some people achieve results more or less quickly than the average number of sessions that we recommend. Long-lasting results cannot be expected after just one or two sessions. Here again, the duration of the resting cycle of the hair follicles varies from one area of the body to the next.

Note: The IPL Hair Removal appliance is not effective on white or light coloured hair.

6.6 Using the IPL Hair Removal safely

Safety is a priority with the IPL Hair Removal

Intense Pulsed Light(IPL) technology provides maximum safety and minimum risk.

Epilation using the Intense Pulsed Light technology lets you achieve long-lasting results with less power intensity compared to other photon-epilation methods. The Hair Removal's low energy use reduces the chances of injury or complications and contributes to your overall safety.

The IPL Hair Removal has a mechanism to protect your eyes.

The IPL Hair Removal is designed with a built-in mechanism(contact sensor) that will protect your eyes. This was designed so that the pulse cannot be emitted if the Hair Removal is aimed at open space.The safety mechanism only allows the appliance to work when the application surface is touching the skin. Do not look at the area being treated when the unit flashes.

7. Use


The IPL Hair Removal offers a safe and long- lasting epilation in the privacy and comfort of your own home.

The IPL Hair Removal is a hair-removal appliance that eliminates unwanted hair from the body(in particular,legs and arms), the sensitive areas (underarms,bikini line) and from the face (cheeks,upper lip and chin).

The IPL Hair Removal can be used by men and women aged 18 and older.

If you are unsure about using the IPL Hair Removal , please consult your dermatologist or your GP.

7.1 Description of the IPL Hair Removal appliance

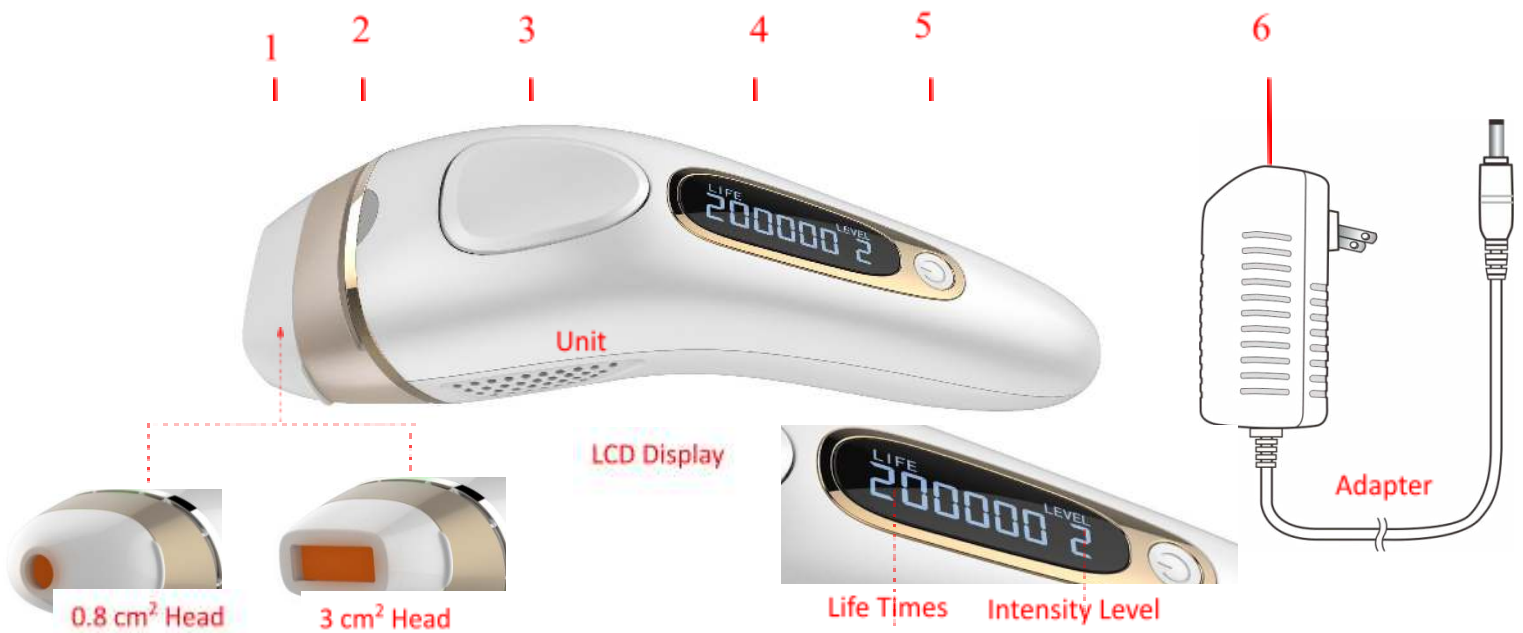

1. Contact sensor heads (0.8cm² & 3cm²Head)
2. Ready indicator
3. Release mechanism
4. LCD Display
5. Button: ON/OFF + power level intensity selector
6. adapter

7.2 Prepare your skin for treatment

Your skin should be shaved, clean, dry and free of powder and deodorant residue.

Do not epilate (whether using an electric Hair Removal or wax) before using the IPL Hair Removal.


The epilation (removal of hair at the root) will make the pulsed treatment ineffective.

7.3 Installing and configuring your IPL Hair Removal appliance



7.3.1 Remove the appliance and other IPL Hair Removal parts from the box.

7.3.2 Plug the DC head of the adapter into the appliance.

7.3.3 Plug the plug of the adapter into the mains.

7.3.4 Press  button to switch the appliance on. The LCD Display and Ready light will light on, and you will hear a slight noise of a fan.

7.3.5 Your IPL Hair Removal is ready to use.

7.3.6 Select the power level that corresponds to your skin type (skin card). Using the  button, select the desired intensity level. To move up from level 1 to level 4, for instance, press the  button as often as necessary, the appliance will move up from level 4 to level 5, then backwards from level 5 to level 4, from level 4 to level 3 and so on down to level 1.

7.3.7 Press the application surface firmly against your skin so that there is full contact. When the application surface is in full contact with the skin. The READY indicator light will start to flash.

NB: The appliance has a mechanism to protect your eyes, a contact sensor that prevents the appliance from releasing a pulse if its application surface is not in full contact with your skin. If the appliance is not in full contact with your skin, it will not release a pulse of light. Make sure, therefore, that application surface is correctly applied to your skin.

7.3.8 Press the release mechanism. The appliance will release a pulse of light, You may feel a slight

warm and tingling sensation.

7.3.9 Remove the application surface from the treated patch of skin.

7.3.10 Move the application surface to another area to be treated. Pressure marks left on your skin by application surface should help you to correctly align the Hair Removal to administer the next pulse, avoiding the release of two pulses on the same patch of skin or too much space being left between the patches (see 7.4.4 'Cover the area being treated') . Your IPL appliance will recharge immediately after each pulse. After 1.6 second (at level I), the READY light will come on again. The appliance is now ready to release another pulse.



Warning: Never treat the same patch of skin more than once per depilation session! Repeated treatment of the same area increases the risk of side effects.

7.4 Modes of use

7.4.1 PRECISION mode

The precision Mode, as described under points 7.3.7 to 7.3.10 above, allows for safe treatment of the accessible and most sensitive areas of your body (face, underarms, bikini line, etc.).

7.4.2 GLISS mode

When the applicator window is properly applied to the skin and the release mechanism is held in, the user can activate the Gliss mode (continuous release of 36 pulses per minute, at level I). In this way, there is an automatic continuous release of pulses without the need for additional handling. You can thus easily slide the appliance over your skin without interruption in a slow and continuous movement without lifting or lowering it.


Note:

If you lift the appliance while holding the release mechanism pressed in, the Gliss mode remains active. But for your safety, no pulse is released until the appliance is in full contact with your skin. Ensure that the contact sensor that forms the contour of the appliance's applicator surface is perfectly and firmly pressed against your skin to guarantee operation of the continuous flashes in the Gliss mode. If you break the contact between the applicator surface and your skin, you have 5 seconds to adjust the appliance to ensure full contact with the skin while holding the release button pressed in. If, after 5 seconds, contact is not full, the Gliss mode will deactivate and will have to let the release button go and press it again to restart the process.

If you lift the appliance and let go of the release mechanism, the Gliss mode will be deactivated.

The Gliss Mode is suitable for people who already have some experience and some dexterity in using the appliance. With 36 pulses per minute (at level I), the Gliss Mode lets you treat larger areas of your body more safely and more quickly: legs, back, arms, torso, ect. For the sensitive areas and the face, use the Precision Mode.

7.4.3 Automatic shut-off feature

If the appliance is left for more than 10 minutes, it will automatically shut itself off, The 'READY' indicator light will then be orange. To use the appliance again, the user must press the  button and reset the appliance to the desired power level.



Avoid releasing more than one pulse on the same patch of skin!
Do not treat the same area more than once in one session!

7.4.4 Cover the area being treated

Pulses must be administered in rows, from one end to the other and proceeding Systematically in successive rows. The method ensures better control of pulses and will help you to cover the area better while preventing treating the same Patch twice or administering pulses too close to each other in the same area. Pressure marks left on your skin by the application surface should help you to correctly align the Hair Removal to administer the next pulse. Avoiding the release of two pulses on the same patch of skin or too much space being left between the patches.



Stop the treatment immediately if your skin has burns or blisters appear!

7.5 Epilation of the face with the IPL Hair Removal.



It is imperative not to look directly at the area treated during the flash and to wear protective eye glasses when using the IPL Hair Removal to treat the face.
Do not use it for epilation of eyelashes, head hair, and eyebrows.
This product is not appropriate for full removal of male beard.



For your comfort and safety, we recommend you have someone help you for the first few sessions using the IPL Hair Removal on your face.
The PRECISION mode is required when removing facial hair.
DO NOT USE THE GLISS MODE TO TREAT THE FACE OR SENSITIVE AREAS.



To treat the face, it is imperative that you wait 5 seconds between flashes. The duration of the IPL Hair Removal depilation session on the face should not exceed 10 minutes.

7.6 Approximate duration of the session.

Area	Average duration	Mode
Half-leg	8 to 10 minutes	Gliss
Leg	15minutes	Gliss
Bikini line	2 minutes	Precision
Underarm	1 minute	Precision

7.7 What to expect with the IPL Hair Removal

Using the IPL Hair Removal at home could for people be a new experience in a depilation method that uses light, The IPL Hair Removal is easy to handle and the depilation sessions are quick. You can expect to see, hear or feel the following during a session.

Fan noise – IPL Hair Removal makes a slight noise of fans, This is normal.

Moderate pressure on the application surface – necessary and useful to position it correctly over the connecting area; this is an exclusive safety feature of the IPL Hair Removal.


Heat and prickling sensations – Whenever you release a pulse, it's normal to feel some heat and prickling resulting from the energy of the light.

Do not forget that it is important to use the lowest power setting during the first session. The heat sensation may persist for up to an hour after your depilation session.

Skin slightly red or pink – Just after use or within 24 hours of using it, it is not rare to see the skin turn slightly red or pink, especially around the hair.

If obvious redness, blisters or burns appear, stop using the IPL Hair Removal immediately.

7.8 After a treatment with your IPL Hair Removal.

- At the end of your session with the IPL Hair Removal, switch your unit off using the  button.
- Unplug the adapter from the mains.
- We recommend cleaning your IPL Hair Removal after each session, particularly the application surface (see 'Cleaning your IPL Hair Removal').
- After cleaning your appliance, we recommend storing it in its box and away from any source of water.



Care of your skin after you depilation session

Avoid exposing the treated areas to the sun. Protect your skin with a sun throughout the treatment and for at least two weeks after the last session with your IPL Hair Removal appliance.
Do not use scented cream or deodorant on the treated areas after your sessions.



Side effects and complications.

Some patients may notice a change in pigmentation following a treatment with the IPL Hair Removal. These effects, if they appear, are generally slight and temporary.

If you experience complications (see 'Possible side effects'), stop using the appliance immediately and contact your GP.(General practitioner).

8. Maintenance of your IPL Hair Removal

8.1 Cleaning your IPL Hair Removal.

We recommend you clean your IPL Hair Removal after every session. Particularly the application surface.

Unplug your appliance before cleaning.

Gently clean the surface of the appliance and particularly the application surface with a clean dry cloth and a cleaning agent formulated for cleaning electronic equipment..

Never immerse your IPL Hair Removal appliance or one of its parts in water.

8.2 Autonomy of your IPL Hair Removal

With an exceptionally long life, the IPL Hair Removal flash lamp releases all its light pulses without a decline in energy.

When the appliance is at the end of its life, the 'READY' indicator light will turn red.



Do not use the unit if large stains appear inside the flash lamp or if the application surface is broken. Contact the customer service.


9. Troubleshooting

9.1 My IPL Hair Removal won't start.

Ensure the adapter is properly connected to your IPL Hair Removal and plugged into the mains.

9.2 There are no pulses when I press the release mechanism.

Ensure that the application surface is evenly and firmly in contact with our skin and that the 'READY' indicator light is flashing ; this indicator light guarantees that the appliance is in full contact with your skin, For your safety, the release switch will be able to release a pulse only if the application surface is firmly pressed against your skin.

- Make sure the 'READY' light is lit.
- If the 'Ready' light is still off after 10 seconds, switch your IPL Hair Removal off and again by pressing the  button twice. If the problem persists, contact the customer service in your area.
- The 'READY' indicator light is red and flashing. The appliance has switched to the safety mode because it has detected that it is overheated or there is a critical problem (broken fan, capacitor problem, etc.). Under circumstances like these, switch the appliance off and allow it to cool for 15 minutes before switching it on again.
If the problem persists, contact the customer service in your area.
- The 'READY' indicator light will flash green or red alternately. The adapter being used to connect your IPL Hair Removal appliance is not supplied with your appliance or it is defective. If the problem persists, contact the customer service in your area.
- When the display screen shows '000000' and 'READY' indicator light is red, the appliance has run out of its useful time.



Do not attempt to open or repair your IPL Hair Removal yourself.

Only certified repair centers are authorised to make repairs.

Opening your IPL Hair Removal could expose you to dangerous electrical parts and the energy of the pulsed light, which could cause serious physical injury and / or irreversible eye injury.

Opening your IPL Hair Removal could also damage the appliance and will void the warranty.

Please contact customer service if your appliance is broken or damaged and requires repair.

10. Customer Service.

If your device is broken, damaged, needs repair or if you need any assistance, please contact customer service.

11. FAQ

11.1 Is the IPL Hair Removal truly effective?

Yes, Clinical tests have proven that the IPL Hair Removal achieves long – lasting results in complete safety.

11.2 On what parts of the body can I use the IPL Hair Removal?

The IPL Hair Removal has been designed for removal of body hair. The areas most commonly treated with the IPL Hair Removal are the legs, underarms, arms, the bikini line and the face (cheeks, upper lip and chin).

11.3 How long is a depilation session with the IPL Hair Removal?

See chapter 6.6, about the approximate duration of the session.

11.4 Is the IPL Hair Removal safe?

Your safety was our first priority when we designed the IPL Hair Removal, Very strict standards were followed when testing and approving it for use at home skin tolerance tested under dermatological control. But like any electronic appliance that comes in contact with the skin, the user must follow the instructions and scrupulously follow the safety precautions.

11.5 Are treatments using the IPL Hair Removal painful?

When used properly, most IPL Hair Removal users report a slight heat sensation when the unit pulses. Users with darker and thicker hair sense a little more discomfort. But these sensations ease and disappear at the end of the session. For your comfort and safety, the IPL Hair Removal has 5 power levels depending on your skin type.

11.6 At what frequency should I use the IPL Hair Removal?

The first 4 sessions with the IPL Hair Removal should be 2 weeks apart. Next sessions should be 4 weeks apart until you achieve the desired results.

11.7 Is the IPL Hair Removal effective on white, grey or blond hair?

The IPL Hair Removal is most effective on dark hair or on those containing more melanin. Melanin, the pigment that gives colour to skin and hair, absorbs optical energy. Black and dark brown hair respond best to the treatment. If brown and light brown hair also responds, these will need a few more sessions. Red hair may also respond partially to the treatment. Generally, white, grey or blond hair does not respond to the treatment, but some users have noticed results after several depilation sessions.

11.8 Can I use the IPL Hair Removal on brown or black skin?

Do not use the appliance on naturally dark skin! The IPL Hair Removal eliminates hair by targeting the follicular pigment. Pigments are also found in varying quantities in the surrounding dermal tissues. The quantity of pigments in a person's dermis, visible through the colour of the skin, defines the level of risk to which he or she is exposed by using the IPL Hair Removal. Treating dark skin with the IPL Hair Removal may involve risks like burns, blisters and changes in skin colour (hyper – or hypopigmentation). Please consult the table showing the different phototypes and the use recommended according to these types in the section 'Use'.

11.9 When will I see the first results with the IPL Hair Removal?

Like all other photo-epilation methods, the results are not immediate; you may even have the impression that nothing is happening. The hair will sometimes appear to continue growing after a

session, but after two weeks, most of this hair will simply fall out. This is why you will need to have several sessions before achieving the expected result.

11.10 Can men use the IPL Hair Removal?

The IPL Hair Removal can be used by men. However, men's hair, especially chest hair will require more sessions to achieve the desired results.

11.11 Why does hair grow back in the area that I treated last week?

It is quite common for hair to appear to continue to grow for one to two weeks after a depilation session with the IPL Hair Removal. This process is known as 'ejection'. After two weeks, you will nonetheless notice that this hair falls out or comes away from its follicle. That said, we recommend you do not pull the hair from the follicle - let it fall out naturally. Moreover, some hair will not be affected by the IPL Hair Removal either because of poor application or because the hair was in its dormant phase. This hair will be treated in the following sessions, hence the need for several sessions to achieve the expected results with the IPL Hair Removal

11.12 I've heard that some hair grows back lighter and finer after photo-epilation. Is this true?

This phenomenon is well documented by beauticians and doctors who use pulsed light appliances for epilation treatments. It is possible that some hair grows back lighter and finer after an IPL Hair Removal session. Generally, this hair only represents a tiny part of the original hair growth, treatments should therefore be continued to achieve the desired result.

11.13 Is long-term use of the IPL Hair Removal dangerous for my skin?

Over the last 15 years, the use of the pulsed light in aesthetic medicine has been well documented in professional journals. These journals and institutions have not reported any side effects or injuries connected with long-term use of pulsed light appliances.

11.14 Can I use the IPL Hair Removal to eliminate chin or facial hair?

IPL Hair Removal allows facial hair removal (cheeks, upper lip, and chin). However, the IPL Hair Removal may not be used for epilation of eyelashes, eyebrows or head hair. Please wear a pair of protective glasses when using.

11.15 How long should I wait before using the IPL Hair Removal if I am not protected from the sun's rays?

You should avoid exposure to the sun before and after treatment with the IPL Hair Removal no matter which part of the body has been treated. You must use suitable protection. e.g. a covering garment or a sun cream (with a protection factor of 50+) which you should then use regularly for three weeks before and after treatment. Ideally, you should wait 4 weeks before using the IPL Hair Removal if you have been exposed to the sun without protection. If you have the least doubt about exposure to the sun consult your doctor or contact the customer service.

11.16 What should I do before I use the IPL Hair Removal?

Before each session with the IPL Hair Removal, it is important that the area to be treated has not been exposed to the sun for at least four weeks. A sunscreen with a high protection factor (protection factor of 50+) could be of some help, as well as clothes that cover the area to be treated. Moreover, the area to be treated should be washed with soap and water beforehand and then shaved clean of hair.

11.17 How do I take care of the treated area after using the IPL Hair Removal?

The area treated with the IPL Hair Removal may be cleaned and cared for with well-known skin care products. Specific precautions should be taken to avoid unprotected exposure to the sun. Sun screens with a higher protection factor (protection factor of 50+) and long clothing are good protection against the sun.

11.18 Should I remove hair before, between and after treatments?

No, do not remove hair before, during or after treatments. The hair will fall out on its own naturally and gradually. This stage may last up to two weeks.

12. Specifications.

SPECIFICATION			
UNIT		POWER ADAPTER	
Light source	Quartz bulb	Input:	AC100-240 50/60Hz
Wavelength	550-1200nm	output:	12V <u> </u> 3A
Input	12V <u> </u> 3A		
Area of impact	Interchange head-1	0.8cm ²	
	Interchange head-2	3.0cm ²	
Repetition rate	Level 1	1.6 Sec.	
	Level 5	3.5 Sec.	
Light intensity	Level 1	2.0J/cm ²	
	Level 2	2.5J/cm ²	
	Level 3	3.0J/cm ²	
	Level 4	3.5J/cm ²	
	Level 5	4.0J/cm ²	

ENVIRONMENT			
TEMPERATURE		HUMIDITY	
OPERATING	0 °C-35°C	OPERATING	10% to 60% RH
STORAGE	-5 °C-40°C	STORAGE	5% to 90% RH

This appliance is not intended to be used near a flammable mixture of anaesthetic and air or with oxygen or nitrous oxide.