Manual

Ai08

Beauty your life

PRINCIPLE

IPL photon hair removal uses a specific wavelength, the light penetrates the epidermis reaches the straight skin layer, absorbs a lot of photons to generate heat, and is selectively absorbed by the melanin in the hair follicle. The temperature of the hair follicle rises rapidly destroying the hair follicle, and the photons absorbed by the surrounding tissue are not affected. Injury, the growth of hair is delayed or even completely stopped in a short period of time, so as to achieve the purpose of permanent hair removal. Similarly, it can restore the original elasticity of the skin, eliminate or reduce wrinkles, and reduce pores. Improve skin texture, complexion and firm skin.

HAIR REMOVAL

Using the penetration of strong pulses, it directly reaches the hair follicles, causing the hair follicles to generate instant high temperature, causing the hair follicles to be damaged, and the hair growth rate is slowed down to achieve the effect of hair removal.

PHOTO REJUVENATION

For abnormal pigmentation or abnormal skin blood vessels, photorejuvenation phototherapy uses IPL (strong pulse light) to improve the elasticity of the skin and make the skin look younger. This is due to the IPL promoting the synthesis of collagen in the skin and making the skin look fully elastic and shiny.

PDT

In the process of PDT (photodynamic therapy) treatment of acne, IPL is used as a light source to activate substances to shrink the sebaceous glands.

ABNORMAL SKIN AND BLOOD VESSELS

As other laser treatment methods, IPL can also effectively treat vascular structure abnormalities and skin vascular lesions such as flush tubes.

1

A LOT OF ACNE AND PIMPLES

During the treatment of IPL, a compound porphyrin produced by acne bacteria will appear. When porphyrin absorbs light with a wavelength of about 400 nm, free radicals will be generated, causing damage to acne bacilli. IPL is also effective in treating acne scars.

ABNORMAL SKIN PIGMENTATION

IPL has a significant therapeutic effect on various types of pigmentation abnormalities such as scars, freckles and age spots. Compared with other laser treatments, IPL treatment produces less scabs, so you can apply makeup or wash your face on the day of treatment.

	Laser	IPL
Wavelength	Single wavelength	Many different wavelengths
Treatment	Treatment of 1-2 skin injuries	Treatment of multiple skin injuries
Recovery Time	About 1 week	No recovery time
Pain	Local or general anesthesia	No anesthesia
Side effect	Pigmentation	No side effects
Effect	Effective	Without repeated treatment

NO PAIN AT FREEZING POINT

Ice sensation and hair removal go hand in hand, and enjoy the new ice sensation for painless hair removal.

9-SPEED ENERGY ADJUSTMENT

Adjust the gear according to the state of the skin and hair to remove hair safely and comfortably.

LCD LIQUID CRYSTAL DISPLAY

You can enjoy the beauty salon-grade hair removal at home with a glance at the remaining light-emitting times through the digital display.

2

PRODUCT DESCRIPTION

multifunctional hair removal device (hereinafter referred to as hair removal device) is an innovative beauty product. The hair removal instrument has functions of hair removal, skin rejuvenation and acne removal, and adopts intelligent pulse light irradiation technology, which integrates multiple applications of optical technology in the field of beauty and health. It is one of the most professional and practical products in this field.

PLEASE READ THE INSTRUCTIONS CAREFULLY BEFORE USE

The design of OSENYUAN hair removal instrument is always based on the customer's comfort and safe use. It is manufactured under the highest quality standards. Before using the product, please read and understand all instructions and precautions.

APPLICABLE PEOPLE

Hair Removal Device is an over-the-counter device suitable for men and women to remove excess hair.

APPLICABLE SKIN COLOR

Black skin is prohibited.

Please use the brown skin from the smallest level, and adjust it to the corresponding field according to your skin's ability to adapt to the pain.

3

SAFETY PRECAUTIONS

CONTRAINDICATIONS

The use is prohibited in the following situations:

Children and elderly people, local or systemic inflammation; people with photosensitive skin and people who are using photosensitive drugs, users of afferent medical electronic devices such as heart pacemakers, artificial heart and lung and other life-sustaining devices may cause electronic medical treatment Incorrect operation of the device.

MATTERS NEEDING ATTENTION

1. Pregnant women, under 11 years old, black moles, tattoos, etc. are not recommended. It is not recommended to use sunbeds or sunbeds and sunless tanning products. If you are not sure whether your skin is tanned or sunburned, we recommend waiting a day before using it.

2. Do not use tattoos, dark brown and black spots, these parts may absorb a lot of light and cause skin discomfort.

3. It is forbidden to use around the ears, neck, nipples, genitals or anus. These parts of the skin are sensitive and vulnerable to damage.

4. Do not stay in the same position for too long, or emit light multiple times in one position, which may generate too much heat and make the skin uncomfortable.

If the skin is red and swollen after use and lasts for more than 24 hours, please suspend the use, please contact the doctor or consult customer service.

SAFETY WARNING

1. There is high voltage in the hair removal instrument, and it is strictly forbidden to disassemble it without permission.

2. When using the epilator, please do not cover the vents.

3. The hair removal device must not get wet or water.

4. It is forbidden to use the epilator when the plug or power cord is damaged.

5. It is forbidden to operate the epilator in flammable and explosive environments.

SAFE OPERATION SPECIFICATIONS


1. During flash output, it is forbidden to point the bald head at the fragile parts such as the eyes. Wear goggles when using.

2. Before using this product, please remove all metal objects carried on the body, such as ring bracelets, etc., so as not to cause equipment failure.


3. Please operate in a place away from children. After use, it should be placed in a place where children can't reach it.

4

BUTTON FUNCTION INTRODUCTION

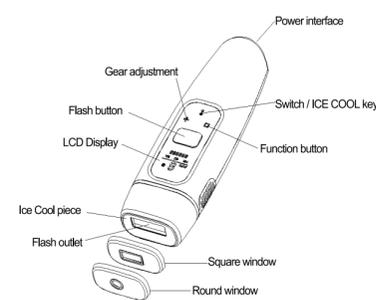DISPLAY INTRODUCTION

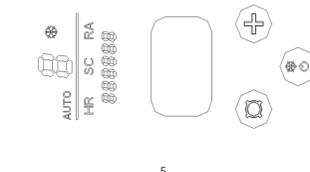5

LCD DISPLAY	FEATURES	PARAMETER
HR	Hair removal	60000 flashes
SC	Skin rejuvenation	
RA	Acne removal	
AUTO	ICE Cool	
☼	Automatic mode	/
1-9	Intensity display	1-9
888888	remaining times	/
☼	Function switch button, it can switch between "HR" "SC" and "RA"	
☼	Long press for 2 seconds to turn on/off, short press to turn on/off IEC	
+	Gear Switch, 1st-9th	

SKIN SENSITIVITY TEST


If you are using this product for the first time, perform skin sensitivity test according to the following requirements:

- Turn on the computer according to the following steps.
- Select the gear 1st gear or 2nd gear.
- Operate the selected skin area according to the operation instructions and precautions.
- After 30 minutes of skin test, check the skin reaction.
- If the skin of the test site is normal, you can continue to use it. If there is abnormally, please stop using it.

6

STEPS FOR USAGE

- Before use, please shave the hair off the hair removal area, clean and wipe it. Make sure that the hair follicles under the skin can absorb the energy to the maximum to achieve the best hair removal effect.
- Connect the power adapter, turn on the power, and wear the exclusive goggles equipped with the product. The backlight of the epilator display is on, and the display status is "OFF".

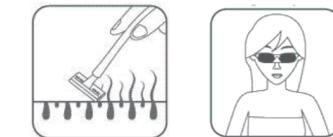- Press and hold the touch button "on/off/ice cool" for 2 seconds, the hair removal instrument starts, refer to the previous "display introduction" to confirm the current working status
- Press the "function switch" touch key, when the screen displays HR, enter the hair removal mode
- The machine is equipped with 1-9 gears. Short press the "+" touch key to adjust the intensity of light wave. For the first time, please use it from the 1st gear. Strengthen the gear according to your own experience. The higher the gear, the better the effect.

7

6. Light setting

Manual light setting (applicable to hair removal in small areas)  
By default, the hair removal instrument starts by manual light emission. Press the flash emission button once to emit light once.

Sliding light setting (applicable to large area hair removal)

After the hair removal device is turned on, press and hold the flash emission button, the hair removal device switches to sliding light, and the "AUTO" on the display will blink continuously. You only need to make the lamp head close to the skin vertically, and the hair removal device will automatically emit light.

7. Close the lamp head to the part where hair removal is needed and use it vertically

(1) When the hair removal device is in manual light output mode, wait for the three touch buttons to flash, press the flash emission button to emit flash, the hair removal area of the lamp cover is completed, the lamp head moves to the next hair removal area, and close to the skin, wait three The touch buttons start to flash, press the flash emission button to emit flash, and cycle in sequence, it is recommended to irradiate each part 2-3 times.

(2) When the hair removal device is in the sliding light mode, long press the flash emission button, a flash will be emitted, and the area covered by the lamp head is depleted. The lamp head moves to the next hair removal area and is close to the skin. It is recommended to irradiate each part 2-3 times.

Tips: In sliding light mode, the lamp head fully adheres to the skin, and it can flash continuously for 8 seconds when the lamp head slides.

8. If you want to use the ice pack function, simply press the "Switch On/Off button", the ☼ will appear on the display, press the "Switch On/Off button" again, the ☼ on the display will disappear. Then turn off the ice pack function.

9. After use, press and hold the "switch on/off button" for 2 seconds to turn off, the epilator is in the power-on state, no action for 5 minutes, and automatically shut down.

10. Unplug the adapter, and then remove the adapter from the power base.

8

11. After ten minutes, after the surface of the lamp cap has cooled down, wipe the light window of the lamp cap with a non-alcoholic wipe or soft paper towel.

12. After cleaning, put the epilator and adapter back into the packaging box, place the epilator in a ventilated and dry place, and keep it out of the reach of children.

13. Care for the skin in the light-receiving area.

20 minutes after the end of hair removal, moisturizing skin care products can be used to repair the hair removal area. If there is discomfort, cold compress for 30-60 minutes until the skin does not feel burning.

NOTE

1. After removing hair, apply aloe vera gel or similar moisturizing cream to the hair removal area. Do not use harsh cosmetics immediately.

2. Be sure to clean the treatment area before hair removal. It is recommended to wash the skin with cold water. Do not touch the water within 2 hours after the operation. Do not sunbathe, warm water bath, steam bath, etc. within 48 hours.

3. Avoid exposing the hair removal site to strong sunlight within 2 hours after hair removal and apply sunscreen before outdoor activities. It is recommended to operate at night.

4. Do not swim in the disinfected swimming pool within 48 hours after hair removal.

5. Do not press the trigger button continuously, just press it once. It is recommended to irradiate each part 2-3 times. (The area of the skin where the entire lamp head is in close contact is one place).

6. The continuous working time of the machine should not exceed 25 minutes. If the hair removal area is relatively large, it can be divided into several times (the power should be turned off after 10 minutes, and you can continue to use it after 10 minutes of rest)

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